



Bencox[®] ID Stem



Bencox[®] ID Stem

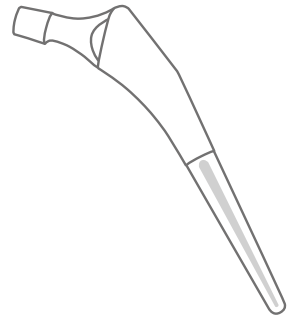

Bencox ID Stem Specifications

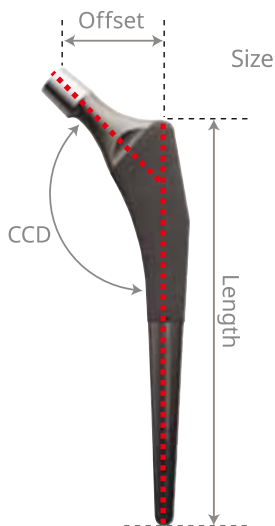

Size	Length (mm)	Offset (mm)	CCD (degree)	Rasp
1	137	37.4		1
2	141	38.1		2
3	145	39.3		3
4	149	40.0		4
5	153	40.7		5
6	157	41.4	132°	6
7	161	42.1		7
8	165	43.3		8
9	169	44.5		9
10	173	45.7		10
11	177	46.9		11

Unified Instrumentation


- Cementless and cemented stem of ID system use same instrumentation.
- This instrumentation maximizes efficiency in the operating room, while still providing intraoperative flexibility.



Proven Ti-Plasma Spray Porous Coating

Clinically proven Ti-Plasma spray porous coating

- Provides the primary and secondary stability.
- Protects osteolysis from the wear debris.

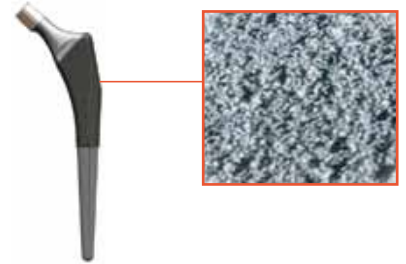

Wide Range of Prosthetic ROM(P-ROM)


Enhanced neck shape

- Patented trapezoidal neck design allows wide angulation and protects impingement & dislocation.
- Same neck design with Bencox Stem is applied.




Compatible Tapered 12/14 Trunnion

- Adoption of Bencox Forte & Delta, and Metal heads
- Provides 3 optional offset (Short, Medium, Long).
- Extra long offset provided for Metal heads

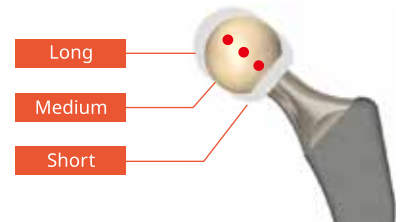


Normalized Proximal Stress

- Proximally tapered stem geometry converts the shear forces to compressive load.
- This helps the stem to resist subsidence and medial migration.

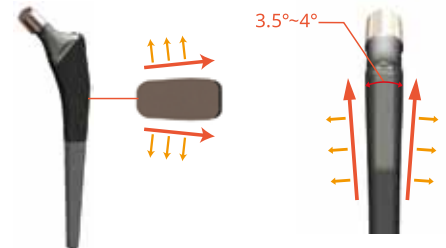

Extensive Stress Distribution

- Low elastic modulus of Ti alloy and the double tapered stem reduce stiffness which avoids distal load transfer and prevents proximal stress shielding.

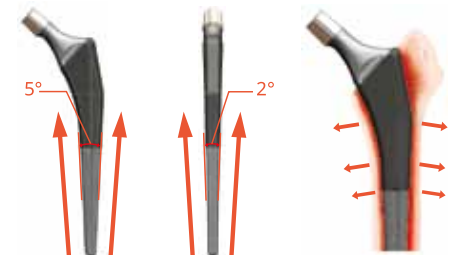

Enhanced Rotational Stability

- Rectangular cross-section
- Longitudinal grooves
 - Longitudinal grooves are designed to resist rotational motion.



Bencox[®] ID Stem ●

Ordering Information

Part No.	Description
01.05.021	Stem 1
01.05.022	Stem 2
01.05.023	Stem 3
01.05.024	Stem 4
01.05.025	Stem 5
01.05.026	Stem 6
01.05.027	Stem 7
01.05.028	Stem 8
01.05.029	Stem 9
01.05.030	Stem 10
01.05.031	Stem 11

Manufacturing Plant 12, Yeongsanhong 1-gil, Ipjang-myeon, Sebuk-gu, Cheonan-si, Chungcheongnam-do, Korea 31056
Tel : 82-41-585-7114 Fax : 82-41-585-7113

Seoul office 8F, 11F, Chungho Tower, 483, Gangnam-daero, Seocho-gu, Seoul, Korea 06541 Tel : 82-2-3445-5492~5 Fax : 82-2-3445-5467,5497

Corentec America 60 Washington Street, Suite 202, Morristown, NJ 07960, USA Tel : 1-862-777-8320 Fax : 1-862-777-8319

EC Representative : Ergo Europe Molenstraat 15, 2513 BH, The Hague, The Netherlands Tel : 31-70-345-8570 Fax : 31-70-346-7299